Micro Inspector


Coliform

Fast · Accurate · Easy To Use Self-filling Test Ampoule

Water Test Kit

Click to return to swimming beach & marine water testing & sanitation information at InspectAPedia.com

Original Source:
Simpltek Water Testing, 8 Byron
Rd. Commack, NY 11725, Tel:
1-631-864-0099, Email:
info@simpltek.com, or
sales@simpltek.com, Website:
http://www.simpltek.com

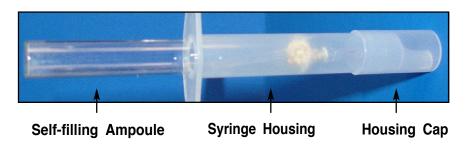

Overview

The Micro Inspector is a unique, self-contained single-use test that screens for aerobic microbial activity. Test results indicate presence/absence and activity levels measured in cfu's. The self-filling glass test ampoule is housed in a protective plastic tube. Ampoule filling is accomplished by simply pushing down, and breaking off ampoule tip in water sample.

Operation Steps

- 1. Remove the housing cap.
- 2. Fill supplied water sample cup (contains neutralizing de-chlorinating solution) or place ampoule tip in water sample.
- 3. Technique #1 push down With ampoule tip placed in water sample, grip syringe housing as shown on right and push down until tip breaks against housing side and ampoule fills. Be sure to keep ampoule tip below water line to completely fill ampoule.
- 4. Technique #2 twist
 With ampoule tip placed in water sample, grip syringe housing with one hand and
 ampoule with other hand. Twist and push down lightly until tip breaks and ampoule
 fills. Be sure to keep ampoule tip below water line to completely fill ampoule.
- When ampoule is completely full, remove from water sample and gently rock ampoule back and forth to totally disolve bio-chemical reagents.
- 6. Incubate at approximately 75F degrees for 24 hours.

Observe Results

Observe color at bottom of test ampoule and match to this chart.

Color	Inspection Time		
	4 Hours	12 Hours	24 Hours
START	Negative	Negative	Negative
YELLOW	Dangerous	Problem	Warning

Presence/Absence: At end of 24 hours observe for color change.

Semi-quantitative: Observe at 4, 12 and 24 hours for color change. The earlier the color change the greater bacteria

amount.

What's in the box: part # ST9300MI

Contains: 1 test ampoule in syringe housing, 1 sample cup with dechlorinating solution, 1 Instruction Sheet.